

# Freedom NXT

CREATIVE THINKING



The NXT is the most rigid, folding, paediatric tilt-in-space manual wheelchair on the market. It is strong and durable, yet compact and lightweight. When the growth spurt hits and the frame won't fit, the NXT adjusts easily for a good fit today and tomorrow.



Yes, you can.®





# Freedom NXT

No one has to choose durability over portability. The NXT incorporates the best design features of both folding, and non folding tilt-in-space mobility. It is lightweight, compact and portable yet strong, durable and dependable. The NXT is also quick and easy to adjust to continue to meet kid's needs as they grow.

## Features




**For kids on the go,** we put our best foot forward...The foot-activated tilt and wheel locks allow easy and secure positioning.



**Transforms as easy as 1, 2, 3...**

- 1: Pop off the wheels and footrests.
- 2: Flip back the armrests and fold down the back.
- 3: Pull up on the strap to fold and stow the chair away for transport.

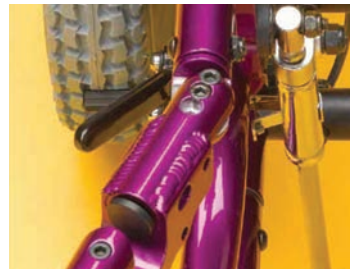

**LEAF (Lower Extremity Accommodation Frame)** is a creative solution to complex seating applications. Unique to Freedom Designs, LEAF is designed so that the front frame angle can adjust to accommodate abduction, adduction and windswept leg positioning.



**MAPS (Multi-Axis Positioning System)** footrest is design innovation for foot and ankle positioning needs. MAPS adjustment feature allows for combinations of dorsiflexion or plantar flexion and inversion or eversion positions, as well as rotation.

## Tech Specs



**Seat width**  
10 - 20"



**Seat depth**  
12 - 24"



**Seat height**  
12.5-18.5"



**Tilt range**  
45°



**Total width**  
seat width +20cm



**Max. user weight**  
100 kg



**Transport weight**  
w/out rear wheels  
10kg



**Max. safe slope**  
6°



**Crash Test Approved**  
Fulfills the requirements of ISO 7176-19

FOR FURTHER INFORMATION:

**New Zealand:** Freephone 0800 468 222, email [sales@invacare.co.nz](mailto:sales@invacare.co.nz) or visit [www.invacare.co.nz](http://www.invacare.co.nz)

**Australia:** Freephone 1800 460 460, email [orders@invacare.com.au](mailto:orders@invacare.com.au) or visit [www.invacare.com.au](http://www.invacare.com.au)

All rights reserved. All information quoted is believed to be correct at time of publication. Invacare reserves the right to alter product specifications without prior consent. 05/2024



**Yes, you can.®**