

Portable Ramps

DEPLOYED IN SECONDS



Invacare Portable Ramps are manufactured with a unique combination of uni-directional glass and carbon fibre materials to create an extremely strong and lightweight ramp suitable for a wide range of uses.



Yes, you can.®

Making Life's
Experiences Possible

Portable Ramps

- **Portable** - Invacare Portable Ramps can be quickly folded into a single compact unit - allowing for easy storage or transport in any vehicle or home environment
- **Lightweight** - Made from a unique combination of uni-directional glass and carbon fibre materials. Invacare Portable Ramps are very strong, yet light enough to be portable
- **Compact** - Multi-paneled design held together with a patented hinge system connected by webbing enabling Invacare Portable Ramps to fold with a concertina action, allowing for a much smaller width when folded
- **Safe** - Clearly marked and grip tape covered, Invacare Portable Ramps are also fitted with pads at the top and base. Invacare Portable Ramps provide a non-slip surface, for all weather conditions
- **Strong** - Invacare Portable Ramps are tested to AS/NZ Standards, up to 300kg!* for security and peace of mind
- **Simple to use** - Invacare Portable Ramps require no assembly allowing the ramp to be deployed in seconds, and immediately ready for use!


Features



Carry bag for Personal and Multipurpose models.



Marked 'TOP' to show head of ramp, grip tape, pads at both ends.

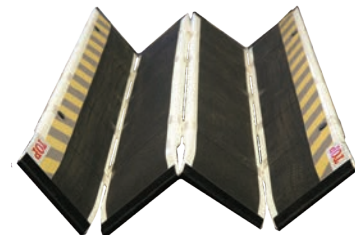

Multi-paneled, hinge system folds with a concertina action, allowing it to be folded to a single compact unit.

Invacare Portable Ramps are available in 4 models, in a wide range of sizes;

- **Personal** - Portable, lightweight ramp that allows easier access in and around the home or workplace. Ideal for climbing kerbs or similar obstacles.
- **Multipurpose** - Longer and wider, designed to cope with a wider range of obstacles making it more suitable for commercial and travel situations. Invaluable when venturing outdoors.
- **Senior** - Ideally used where the doorway is narrow and a great solution for loading unoccupied equipment.
- **Edge Barrier Limiter (EBL)** - A ramp for those who like to get about choose this model for a Ramp developed to give users access to shops, offices and government buildings. It has a 5cm Edge Barrier design for safety & security to prevent steering off the edge of the ramp. This is especially important on the longer length ramps.

Tech Specs

Model/Length cm	Width	Weight
Personal 070	74.5cm	4.7kg
Multipurpose 087	83.5cm	5.6kg
Senior 090	75cm	4.8kg
Senior 120	75cm	8.0kg
Senior 165	75cm	10.8kg
Senior 200	75cm	11.7kg
Senior 300	75cm	21.1kg
Senior 350	75cm	24.7kg

Model/Length cm	Width	Weight
EBL 090	82cm	6kg
EBL 120	82cm	9.7kg
EBL 165	82cm	13kg
EBL 200	82cm	13.7kg
EBL 300	82cm	27.3kg
EBL 350	82cm	31.5kg

FOR FURTHER INFORMATION:

New Zealand: Freephone 0800 468 222, email sales@invacare.co.nz or visit www.invacare.co.nz

Australia: Freephone 1800 460 460, email orders@invacare.com.au or visit www.invacare.com.au

All rights reserved. All information quoted is believed to be correct at time of publication. Invacare reserves the right to alter product specifications without prior consent. 04/2024



Yes, you can.®