

# Softform<sup>®</sup> Premier Active 2

WITH RX CONTROL UNIT



A DYNAMIC THERAPEUTIC  
SUPPORT SURFACE



Yes, you can.<sup>®</sup>

Making Life's  
Experiences Possible



# A dynamic support surface with enhanced pressure relief

With the aim of providing healthcare professionals a dynamic therapeutic support surface - offering additional levels of pressure relief - we have added the RX Control Unit to our leading Softform Premier foam mattress to create the Softform Active 2 RX.

The addition of the RX control unit allows users to deliver optimal client support at the touch of a button via active alternating air cells.

Working together with the clinically proven interface of the Softform Premier foam mattress, these alternating air cells assist in the prevention of pressure ulcers, allowing adequate support.

## **'Very High Risk' Static Surface**

A specially cut, high specification foam interface. Uses effective contouring to assist in microclimate management and pressure redistribution, offering comfortable and effective pressure redistribution for clients at 'Very High Risk' of developing pressure ulcers. The static surface can be easily 'stepped up' with the addition of the new control unit, further assisting in the treatment of severe pressure ulcers.

# A complete patient centric solution



The **Softform Active 2 RX** combines the pressure redistributing qualities of specially cut foam castellations at the patient interface, with air cells beneath which can be inflated as required to provide active therapy.

This construction ensures the patient is always supported by enveloping supportive high specification foam\*, whether the control unit (powering the air cells) is switched on or off. The comfortable high specification foam provides effective pressure redistribution and ensures comfort is maintained at all times.

Softform Premier Active 2 is available in **single** or **king single** to provides a safe, dignified support for patients in need of a wider mattress.

\*(NICE Clinical Guideline 179, September 2014)

# 5

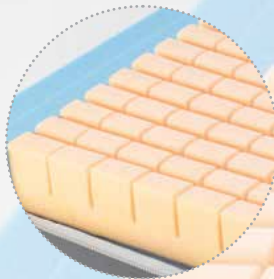
## KEY FEATURES AT A GLANCE

From the Modular Construction to the Plug'n'Play control unit...

### 1

#### High Specification Static Foam Surface

- ▶ The patient will remain on a **'Very High Risk'** pressure reducing surface, minimising the risk of pressure damage in the unlikely event of a power failure.



### 2

#### Plug'n'Play Digital Control Unit

- ▶ An intuitive user interface makes it easy to set up and use.





### 3 Cost Saving Modular Construction

► All components are easily and economically replaceable, **contributing to a lower cost of ownership.**

### 4 Reduce Shear Forces With Glide Mechanism

► An innovative glide mechanism significantly reduces shear forces when used with an electric profiling bed. EPUAP suggests the avoidance of subjecting the skin to both pressure and shear forces, both major contributory factors in the development of pressure ulcers (EPUAP, 2014).

### 5 Self Adjusting System

► No pressure adjustments are required, the RX control unit delivers the optimum pressure.



**Yes, you can.®**

# Powered hybrid support surfaces

There are generally two types of powered hybrid support surfaces.

1

▶ Low density foam within an alternating air cell layer

2

▶ High specification foam layer over an alternating air cell layer

Hybrid support surfaces in 'static' mode ▼



▶ High specification foam layer within a deflated air cell



Hybrid support surfaces in 'alternating' mode ▼



▶ Deflated air cell beneath a high specification foam layer



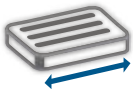
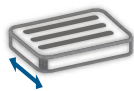
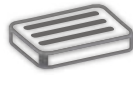


▶ The first type of hybrid support surface is similar to a non-powered hybrid, but offers the additional option of attaching a control unit.





**This control unit has two important functions:**

1. To prevent air from diffusing from the support surface
2. To provide a form of alternating pressure support as a step-up therapy

The alternating pressure support that is being provided has a lower amplitude. This is because the air cells are filled with a low-density foam in the core. Therefore, offloading as an essential part of alternating pressure support is less optimal, compared to traditional alternating pressure support products.

## Tech Specs

	 Length	 Width	 Depth	 Unit weight	 User weight
<b>Softform Active 2 RX Mattress Single</b>	1970mm	880mm	152 mm	14kg	247.6kg
<b>Softform Active 2 RX Mattress King Single</b>	1970mm	1050mm	152 mm	14kg	247.6kg

	 Height	 Width	 Depth	 Unit weight
<b>Control Unit</b>	155mm	275mm	105 mm	1.75kg

## Colour



**Baltic Blue**

**Please note:** colours may vary slightly from those displayed above.

## Other Info

- ▶ High specification castellated foam surface
- ▶ Active support when required
- ▶ Clinically proven for all categories of pressure ulcers
- ▶ Suitable for very high/high risk patients
- ▶ Maximum user weight 247.6 kg

Essential nursing care is pivotal in pressure ulcer prevention. This mattress will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgment and action based planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. A range of assessment scales can be used as a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgment).

For more comprehensive pre-sales information about this product, including the product's user manual, please see your local Invacare website.



**Yes, you can.®**

FOR FURTHER INFORMATION:

**New Zealand:** Freephone 0800 468 222, email [sales@invacare.co.nz](mailto:sales@invacare.co.nz) or visit [www.invacare.co.nz](http://www.invacare.co.nz)

**Australia:** Freephone 1800 460 460, email [orders@invacare.com.au](mailto:orders@invacare.com.au) or visit [www.invacare.com.au](http://www.invacare.com.au)

All rights reserved. All information quoted is believed to be correct at time of publication.  
Invacare reserves the right to alter product specifications without prior consent. **07/2024**



**Yes, you can.®**