

PSP Cushion

POSTURE SEAT POLYMER CUSHION



The Invacare® Matrix® PSP (Posture Seat Polymer) Cushion combines the stability and support of the PS foam base with a polymer overlay for extra protection and shear reduction. An ultra-lightweight cushion providing worry free skin protection and positioning.



matrix[®]
SEATING SERIES

PSP Cushion



Outer cover:

- Reversible Startex™ outer cover can be used smooth or fabric side out
- Moisture resistant and breathable to keep cool and dry
- Includes hook and loop strips both sides

Inner cover:

- Inner cover for extra moisture protection

Foam:

- Molded high resilient polyurethane foam made with Ultra-fresh™ anti-microbial
- Rear cushion radius and waterfall front edge
- Molded polymer insert for extra protection and shear reduction
- Unique waffled relief zone allows immersion in ischial area designed to provide off-loading to hips and thighs for support and stability

Options



Cushion and Cover Modifications:

Standard modifications include seat width and depth reductions, as well as leg length discrepancies



Pelvic Obliquity Kit:

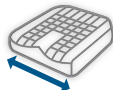

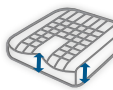



Designed to accommodate up to 1" obliquities.



Cushion Rigidisers:

Lightweight plastic seat insert designed to fit to the bottom of the cushion inside the cover. Sits between seat rails for level support and eliminates hammock effect of sling seat upholstery.

Tech Specs

	 Width*	 Depth**	 Height Adductor/Abductor	 Product weight	 Weight Capacity	 Limited Warranty
Paediatric	14"	12, 14, 16"	3.5" Adductor 3.5" Abductor 2.75" Leg Trough	From 900g (2lb)	68kg (150lb)	2 years
Standard	14, 16, 17, 18, 20"	16", 18", 20", 22"	3.5" Adductor 3.5" Abductor 2.75" Leg Trough	from 1kg (2.2lb)	136kg (300lb)	2 years

* Actual width is 1/2" narrower than stated width
** Actual depth is 1 1/2" longer than stated depth

FOR FURTHER INFORMATION:

New Zealand: Freephone 0800 468 222, email sales@invacare.co.nz or visit www.invacare.co.nz

Australia: Freephone 1800 460 460, email orders@invacare.com.au or visit www.invacare.com.au

All rights reserved. All information quoted is believed to be correct at time of publication. Invacare reserves the right to alter product specifications without prior consent. 05/2024



Yes, you can.®