

Vi Cushion

STABILITY, SUPPORT AND COMFORT



The Invacare® Matrix® Vi cushion is designed to provide skin protection and positioning for individuals at moderate to high risk of skin breakdown. The unique three part foam construction targets skin protection where it is needed most, and features excellent stability, support and comfort.



Vi Cushion



Top Layer:

- Soft High Resilience foam for optimal comfort and support
- Pre-ischial contour provides pelvic stability and abductor/adductor contouring promotes midline leg positioning.
- Made with Ultra-Fresh(TM) for antibacterial and odour protection

Visco foam insert:

- Ischial relief insert made from moulded visco-elastic foam
- Provides immersion to sensitive bony prominences

Bottom Layer:

- Firm High Resilience foam for support and durability
- Trochanteric shelf provides stability to hips
- Waterfall front edge and rear radius for sacral support

Cover:

- Moisture resistant and breathable outer cover including slip resistant base, hook and look fasteners and lifting strap
- Kid*ab*ra print standard in paediatric sizes
- Inner cover is moisture resistant and easily wiped clean

Options



Cushion and Cover Modifications:
Standard modifications include seat width and depth reductions, as well as leg length discrepancies

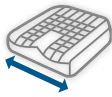
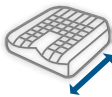
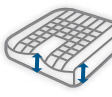





Pelvic Obliquity Kit:
Designed to accommodate up to 1" obliquities.



Cushion Rigidisers:
Lightweight plastic seat insert designed to fit to the bottom of the cushion inside the cover. Sits between seat rails for level support and eliminates hammock effect of sling seat upholstery.

Tech Specs

	 Width*	 Depth**	 Height	 Product weight	 Weight Capacity	 Limited Warranty
Paediatric	12" - 14"	12" - 14"	2.85" Adductor 2.75" Abductor 2" Leg Trough	From 680g (1.5lb)	68kg (150lb)	2 years
Standard	14", 16", 17", 18", 20", 22"	15" - 22"	3.75" Adductor 3.25" Abductor 2.5" Leg Trough	from 1.3kg (3lb)	136kg (300lb)	2 years
Heavy Duty	20 - 30"	18" - 26"	3.75" Adductor 3.25" Abductor 2.5" Leg Trough	from 1.8kg (4lb)	272kg (600lb)	2 years

* Actual width is 1/2" narrower than stated width
** Actual depth is 1 1/2" longer than stated depth



FOR FURTHER INFORMATION:

New Zealand: Freephone 0800 468 222, email sales@invacare.co.nz or visit www.invacare.co.nz

Australia: Freephone 1800 460 460, email orders@invacare.com.au or visit www.invacare.com.au

Asia: Phone +61 2 8839 5330, email anzexportsales@invacare.com or visit www.invacare.com

Thailand: Phone +66 0 2108 8555 or visit www.invacare.com



Yes, you can.®