

Libra Cushion

THE COMPLETE SOLUTION



The Invacare® Matrix® Libra Cushion is designed with the highest level of skin protection, positioning and adjustability for unique individual needs. Specially contoured to provide balance and stability to the pelvis, optimising body function. Recommended for all users including those at highest risk of skin breakdown.



Libra Cushion



Fluid Sac:

- Dual layer, sealed compartments for maximum protection around bony prominences.
- Non-temperature sensitive, no kneading or maintenance required.

Foam base:

- High resilience foam made with Ultra-Fresh™ anti-microbial and odour protection.
- Ab/adductor contour promotes comfortable midline leg positioning
- Pre-ischial contour prevents sliding to reduce shear
- Trochanteric shelf promotes level pelvis and lateral stability
- Waterfall front edge and extended rear cushion radius

Cover:


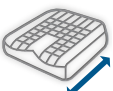
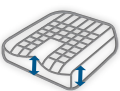



- Moisture resistant, high stretch, breathable Startex™. Hook and loop fastener strips, lifting strap and convenient zippered side pocket.
- Inner platilon liner provides extra lay of moisture protection.

Fit Kit Accessories:

Designed to enhance immersion, correct or accommodate postural issues.

Pelvic Obliquity kit and cushion modifications also available.

Tech Specs

	 Width*	 Depth**	 Height	 Product weight	 Weight Capacity	 Limited Warranty
Libra	14 - 20"	12 - 22"	3.35" Adductor/ Abductor 2.35" Leg Trough	from 1.8kg (4lb)	136kg (300lb)	3 years
Libra Heavy Duty	21 - 24"	12 - 22"	3.6" Abductor/ Adductor 2.6" Leg Trough	from 3.6kg (8lb)	226kg (500lb)	3 years

* Actual width is 1/2" narrower than stated width
** Actual depth is 1 1/2" longer than stated depth

FOR FURTHER INFORMATION:

New Zealand: Freephone 0800 468 222, email sales@invacare.co.nz or visit www.invacare.co.nz

Australia: Freephone 1800 460 460, email orders@invacare.com.au or visit www.invacare.com.au

All rights reserved. All information quoted is believed to be correct at time of publication. Invacare reserves the right to alter product specifications without prior consent. 05/2024



Yes, you can.®