

CS8

FLOORLINE BED



PERFECTLY STYLED IN SINGLE
AND KING SINGLE



Yes, you can.®

Making Life's
Experiences Possible

CS8

EXCLUSIVE FEATURES

One button chair positioning.

Four sectioned electronically operated mattress platform with independent adjustment of backrest angle (maximum 65°) and knee section angle (maximum 35°)

75mm twin Tente castors (four braking castors and four free-wheeling).





Four soft contoured mattress retainer handles.

Removable timber head and foot panels.

Trendelenberg head and foot down tilt 14°.

Under bed lighting.

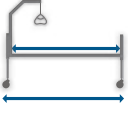

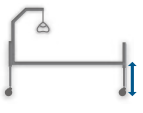


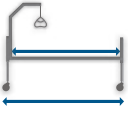

Electrically adjusted height range. The CS8 works within a height range of 150mm low height to a maximum high height of 680mm.



Yes, you can.®

CS8

TECH SPECS

CS8 floorline King Single			
	Length	Width	Adjustable Height
	211cm	109cm	15 - 68cm
CS8 floorline King Single			
	User weight	Bed weight	
	250kg	125kg	
CS8 floorline King Single			
	Mattress Length	Mattress Width	
	197cm	107cm	

Ease of use, functionality, serviceability and cleaning practicality are key features of the CS8.

Basic Dimensions:

Overall Length	211.0cm
Overall Length (base extension extended)	230.5cm
Overall Width King Single	109cm
Overall Width (rails fitted) King Single	120cm
Mattress Surface Length (between panels)	199.5cm
Bed Platform Height Range	15.0cm to 68.0cm
Back Rest Angle Adjustment	0° to 65°
Knee Bend Angle Adjustment	0° to 35°
Trendelenberg Tilt Adjustment (head and foot down)	0° to 14°

Electric Parameter:

Power Input	100-240V AC 50/60Hz, Max current 2.5A
Rate Load	interruption 10%, Max 2 minutes / 18 minutes
Leakage Protection	Class 1 Type B
Waterproof Grade	IP54

Mattress Size:

Recommended Mattress Size King Single	Width = 107cm
	Length = 197cm
Maximum Mattress Depth (when side rails in use)	17.0cm
Minimum Mattress Depth (when side rails in use)	10.0cm

Bed Weight & Safe Working Load:

Without Side Rails King Single	125 Kg (approx)
With Side Rails King Single	144 Kg (approx)
S.W.L Safety Factor = 2 x S.W.L	250Kg

Standard:

CS8 beds are designed and manufactured in compliance with International Standard IEC 60601-2-52 where applicable.

FOR FURTHER INFORMATION:

New Zealand: Freephone 0800 468 222, email sales@invacare.co.nz or visit www.invacare.co.nz

Australia: Freephone 1800 460 460, email orders@invacare.com.au or visit www.invacare.com.au

All rights reserved. All information quoted is believed to be correct at time of publication. Invacare reserves the right to alter product specifications without prior consent. 07/2024



Yes, you can.®